


# Rubber Fab's Tank Vent and Breaker System

# V<sup>2</sup>B GASKETS



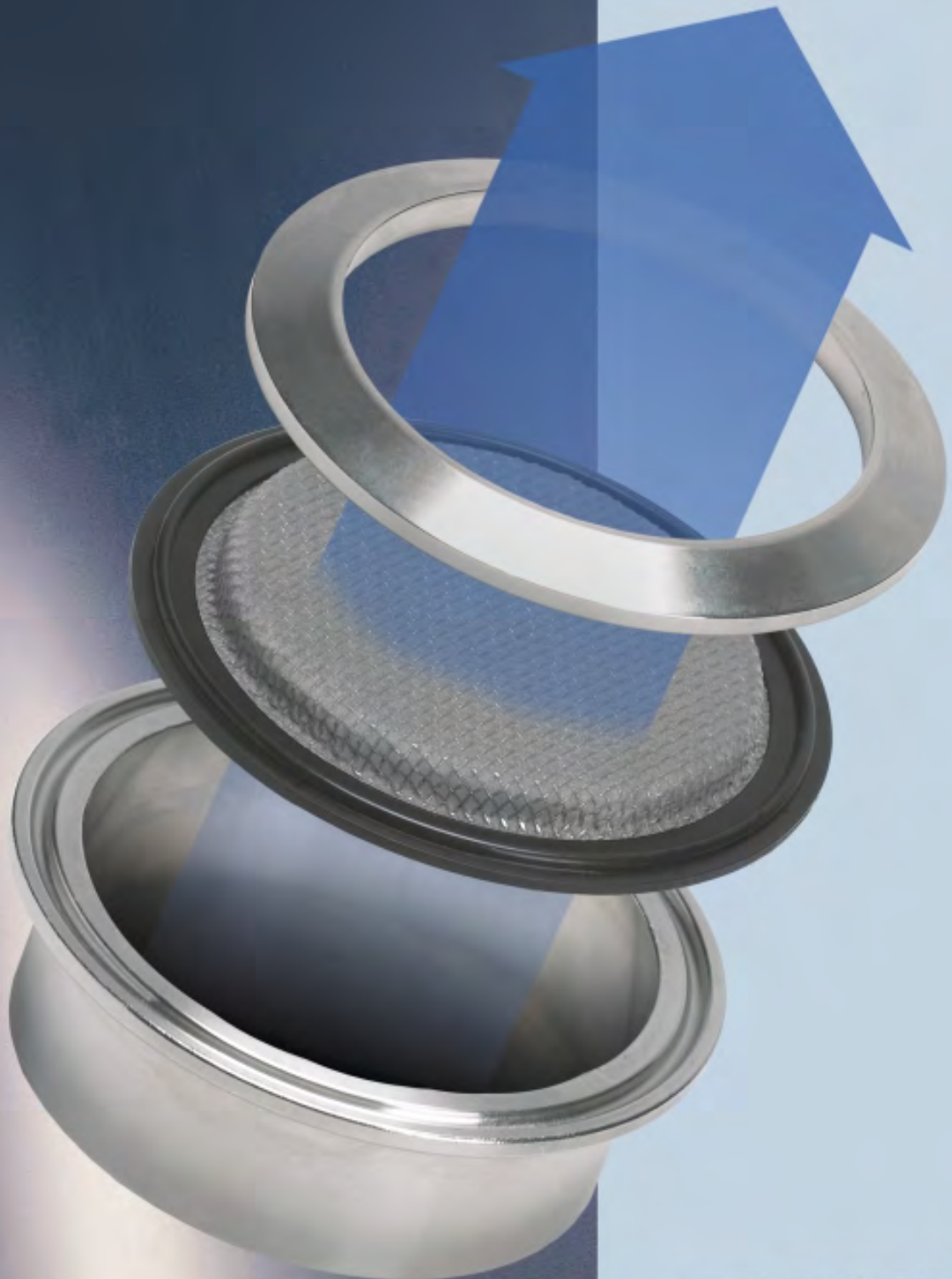

## V<sup>2</sup>B Gasket

- Sizes Available: 1" – 6"
- Removes atmospheric and environmental particulates
- Easy to maintain
- Easy to install
- Meets 3A design
- 316L stainless steel construction
- V<sup>2</sup>B works as both a vacuum breaker and venting system
- Materials of construction are both lot and batch certifiable
- Elastomers available: Platinum silicone, FKM fluoroelastomers, and EPDM
- All V<sup>2</sup>B gaskets are removable and replaceable
- Minimizes need for pre-air filters

## Meets High Pharmaceutical Standards

- U.S. Pharmacopoeia Class VI Certification
- Cytotoxicity Criteria
- Title 21 CFR 177.2600
- USDA Sanitary Standards
- Good Manufacturing Practices (GMP)

## Presenting Rubber Fab's V<sup>2</sup>B Gasket: The Vent and Vacuum Breaker Gasket.

The Rubber Fab V<sup>2</sup>B Gasket offers an additional answer for vacuum and pressure relief for pharmaceutical and food/beverage processing tank venting when airborne particulate removal is required. Using a 200 mesh screen combined with a 14 mesh screen of scintered construction, Rubber Fab has created a cost effective alternative to expensive vacuum breakers and venting systems.

The V<sup>2</sup>B System is attached/connected using a standard sanitary clamp to a sanitary tank port, the same way a standard vacuum breaker is installed. Installation can be horizontal, vertical or in single or multiple port configurations (i.e. manifolds) without special

fabrication or expensive tooling. The V<sup>2</sup>B screen will fit perfectly into the ID of a sanitary pipe and can be used for inflow fluid conditioning filtration.


## Providing Air Flow Without Contamination.

A key benefit of the V<sup>2</sup>B Gasket is it's prevention from most atmospheric and environmental contamination customarily found with other venting systems. Because the V<sup>2</sup>B Gasket can provide unrestricted air flow with ultra fine filtration, you can feel secure in knowing the hygienic condition of the tank contents remain protected. With regular cleaning and maintenance, the V<sup>2</sup>B Gasket is a cost effective and incredibly simple way to provide contaminate free air flow to any tank system.



## Create Greater Air Flow with the V<sup>2</sup>B Gasket.

The Rubber Fab V<sup>2</sup>B Gasket creates the opportunity to increase the flow of air (or decrease if necessary). By using a concentric reducer as illustrated below, you can provide a 2 inch pipe flange 4 inches of air flow area. And by utilizing one or multiple 2-way and/or 4-way sanitary pipe adapters, you can exponentially increase the air flow in and out of your tank system while still protecting the system from atmospheric debris.

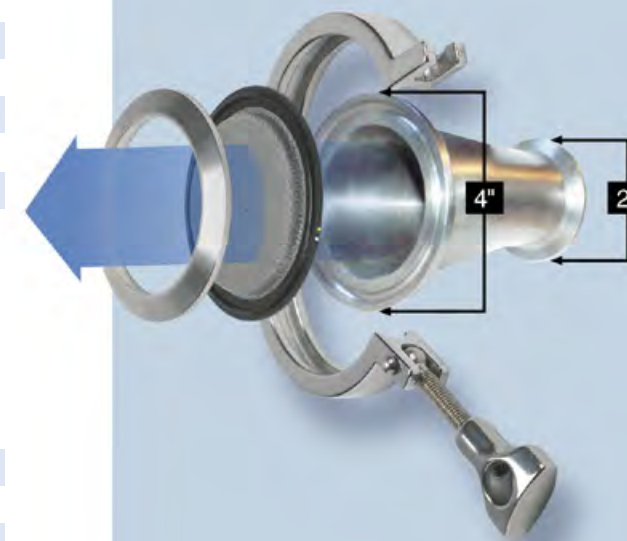

### V<sup>2</sup>B Gasket Flow Rate: Air

Line Size	CFM per unit		
1.5"	1.26 x 6.45* =	8.13 x .51** =	4.15
2"	2.46 x 6.45 =	15.87 x .51 =	8.09
2.5"	4.04 x 6.45 =	26.06 x .51 =	13.29
3"	6.02 x 6.45 =	38.83 x .51 =	19.80
4"	10.92 x 6.45 =	70.83 x .51 =	36.12
6"	25.34 x 6.45 =	163.44 x .51 =	83.35

\* Flow rate through 200 mesh .0021 wire diameter is 930 CFM. When divided by 144, equals 6.45 CFM per square inch.  
 \*\*% open area of 14 mesh .020 wire diameter.

### V<sup>2</sup>B Gasket Flow Rate: Water

Line Size	GPM per unit		
1.5"	1.26 x 16* =	20.16 x .51** =	10.28
2"	2.46 x 16 =	39.36 x .51 =	20.07
2.5"	4.04 x 16 =	64.64 x .51 =	32.97
3"	6.02 x 16 =	96.48 x .51 =	49.20
4"	10.92 x 16 =	174.88 x .51 =	89.19
6"	25.34 x 16 =	405.44 x .51 =	206.77

\* Flow rate of water through 200 mesh .0021 wire diameter is 16 GPM per square inch.  
 \*\*% open area of 14 mesh .020 wire diameter.

### V<sup>2</sup>B Gasket: Open Area

Line Size	Open Area sq/in		
1.5"	1.27 x 1.27 =	1.61 x .7854* =	1.26
2"	1.77 x 1.77 =	3.13 x .7854 =	2.46
2.5"	2.27 x 2.27 =	5.15 x .7854 =	4.04
3"	2.77 x 2.77 =	7.67 x .7854 =	6.02
4"	3.73 x 3.73 =	13.91 x .7854 =	10.92
6"	5.68 x 5.68 =	32.26 x .7854 =	25.34

Please consult factory for CFM/GPM calculations on additional sizes.



The V<sup>2</sup>B Gasket is available in a size range from 1.5" to 6"



## V<sup>2</sup>B Gasket


- Sizes Available: 1-1/2" - 6"
- Available as a **complete system** or as **individual components**
  - V<sup>2</sup>B Screen
  - V<sup>2</sup>B Gasket
  - V<sup>2</sup>B Vent Cap
  - Sanitary Clamp
  - Optional Torque-Rite®
- Use the V<sup>2</sup>B Gasket for:
  - **air flow and venting**
  - **vacuum relief**
  - **particulate removal**
  - **in-line fluid conditioning**
- Used in:
  - **pharmaceutical tanks**
  - **food and beverage processing**
  - **environmental incubation**
  - **clean room applications**
- Components available in electro polished finish.
- **Other mesh and micron sizes are available.** Consult Rubber Fab with special needs and requirements.

NOTE: The V<sup>2</sup>B Gasket is not intended for use in every existing vacuum breaker and venting situation. Rubber Fab Technologies Group recommends testing to determine the individual system's needs and requirements and its compatibility with the V<sup>2</sup>B Gasket. Rubber Fab Technologies Group suggests a regular cleaning cycle in order to maintain consistent water and air flow through the V<sup>2</sup>B Gasket. Due to the inherent risks in vertical silo and tank venting, Rubber Fab Technologies Group does not recommend the V<sup>2</sup>B Gasket without the use of traditional vacuum breakers.

## The V<sup>2</sup>B Gasket: An Excellent Addition to Rubber Fab's Fluid Conditioning Program



Additional information is available on fluid conditioning technology in Rubber Fab's Screen Hygienic Seal brochure. Contact Rubber Fab Technologies Group at 973-579-2959 to receive copies of our literature.

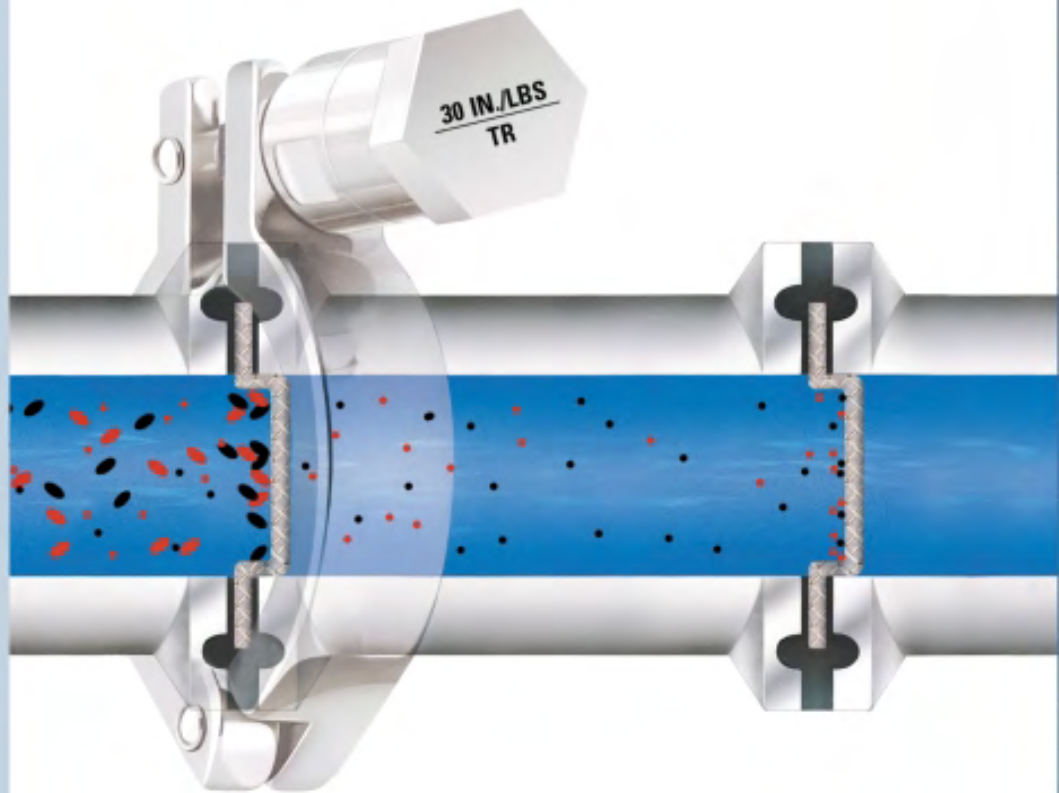

Rubber Fab Technologies Group  
Limecrest Road

P.O. Box 626 • Andover, NJ 07821  
phone: 973-579-2959 fax: 973-579-7275

[www.rubberfab.com](http://www.rubberfab.com)

Rubber Fab Mold & Gasket is a member of:



Distributed by: